


Steriking® Bonded Wrap

For High-Performing Sterile Processing Units.

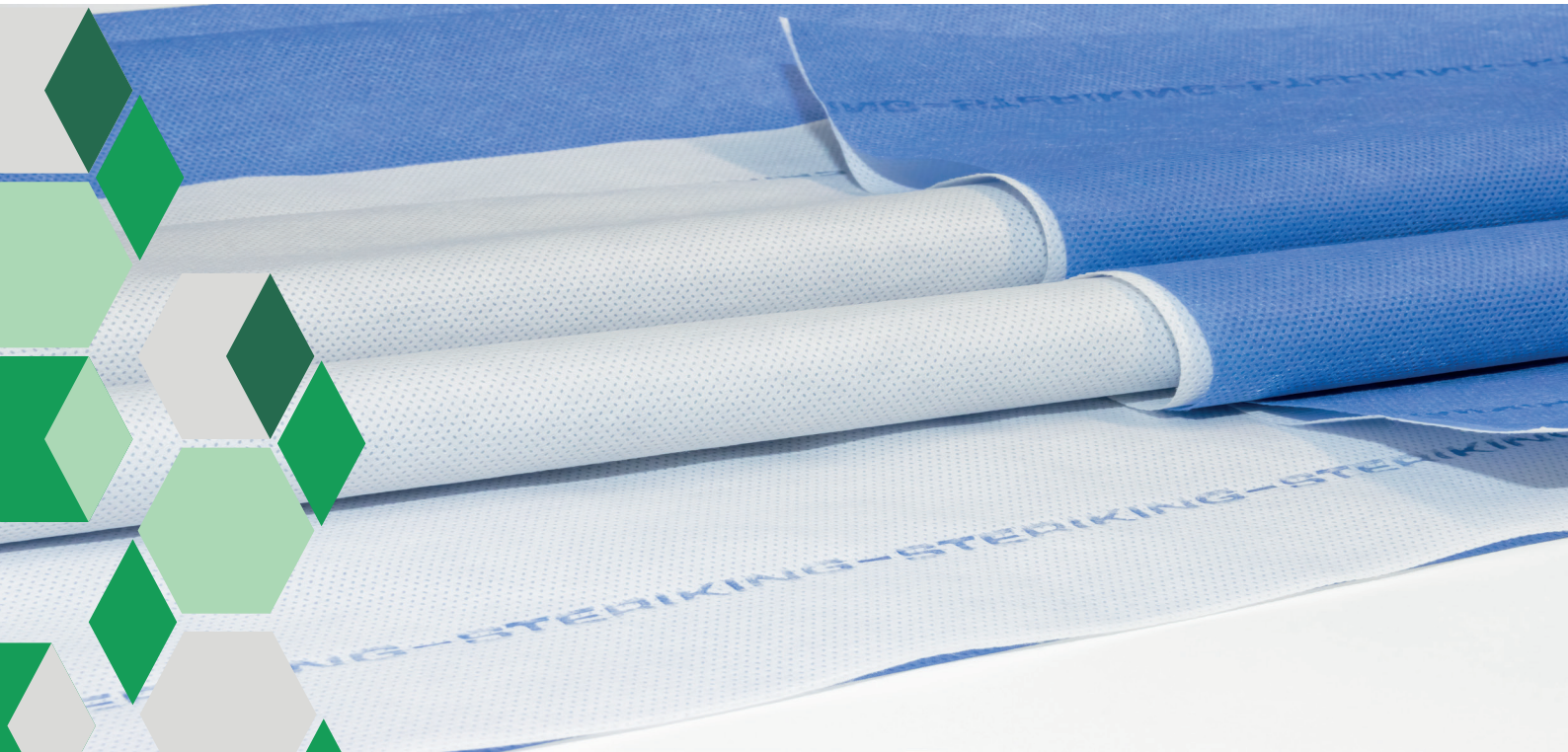


Steriking® Bonded Wraps provide the comfort of easy handling combined with unmatched strength, enhancing operational efficiency in healthcare settings.

Steriking® Bonded Wrap




Steriking® Bonded Wrap products, are crafted from blue and white wrapping sheets which are bonded together. Designed to cater to diverse needs in healthcare facilities, Steriking® Bonded Wrap range offers various material combinations, thicknesses, and sizes for packing instruments, sets, or instrument trays. The cost of packaging depends how it performs. With Steriking® Bonded Wrap, you can ensure patient safety, minimize costs, and minimize the risk of infections and extended hospital stays.

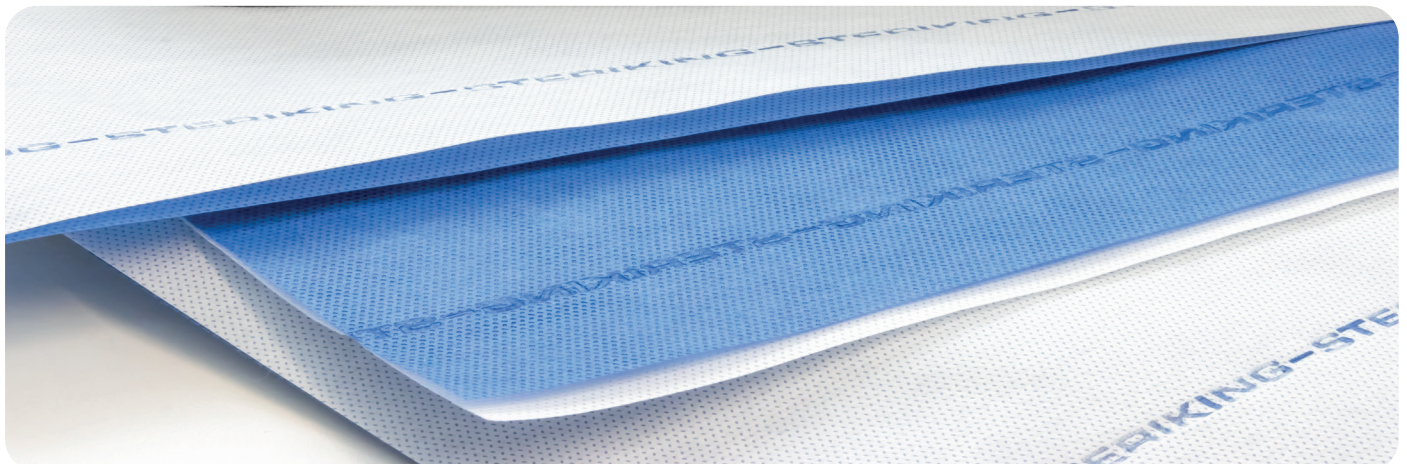

Exceptional Durability and Protection

- Five-layer structure (SSMMS) offers superior resistance against punctures, tears, and abrasions.
- Antistatic properties shield users from static shocks.
- Low linting material prevents particulate contamination.

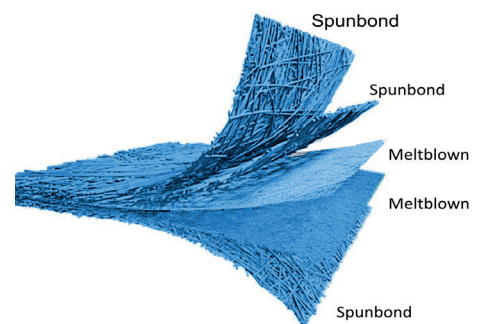

Optimal Softness and Drapability

- Designed with user convenience in mind, Steriking® BW wrap is soft on fingers, with very little folding memory. It's easy to drape, fold, and handle.

Compatible for All Sterilization Methods

- Validated for STEAM, EO, FO, and VH2O2 sterilization processes, ensuring hassle-free selection and reducing the risk of material compatibility errors.
- Reduces the number of needed material combinations.

Material	STEAM	EO	FO	H2O2
Steriking® BW	x	x	x	x

Steriking® Bonded Wrap




Efficient Sterilization, Superior Microbial Barrier

- Bonded wrap material compositions are rigorously tested for sterilant penetration, microbial barrier properties, and bacterial filtration efficiency (BFE) according to industry standards. Ensures a reliable microbial barrier in even the most demanding scenarios.


Colours for Safety and Convenience

- The blue and white bonded sheets, sealed with ultrasonic sealing technique, allow easy identification of any tears, holes, or damage in the wrap.
- The white colour highlights potential impurities effectively.

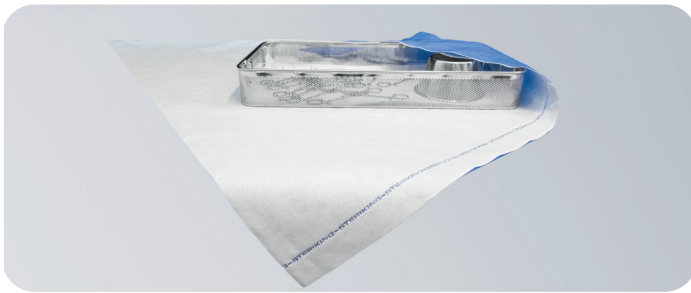

Efficiency at Its Core

- Bonded Wrap reduces the working time in CSSD and saves time during aseptic presentation.
- Wrapping with Steriking® Bonded Wrap takes only half the time compared to the traditional two-layer method.
- Opening of the pack takes only half of time compared to two-layer method.

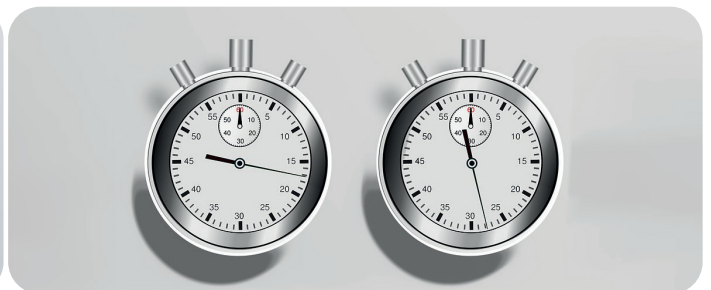

Recycle-Ready

- Made from 100% pure polypropylene, our material supports recycling in plastic streams, contributing to sustainable practices.



Steriking® Bonded Wrap



Steriking® Bonded Wrap BW1- available range

STERIKING® BONDED WRAP 1 (37+47gsm)	
BW1060	60 x 60
BW1075	75 x 75
BW1090	90 x 90
BW1100	100 x 100
BW1120	120 x 120
BW1140	140 X 140

Steriking® Bonded Wrap BW2- available range

STERIKING® BONDED WRAP 2 (47+47gsm)	
BW2060	60 x 60
BW2075	75 x 75
BW2090	90 x 90
BW2100	100 x 100
BW2120	120 x 120
BW2140	140 x 140

Steriking® Bonded Wrap BW4- available range

STERIKING® BONDED WRAP 4 (75+47gsm)	
BW4060	60 x 60
BW4075	75 x 75
BW4090	90 x 90
BW4100	100 x 100
BW4120	120 x 120
BW4140	140 x 140



Wipak Oy
Wipaktie 2
FI-15560 Nastola



We care that you pack safely!
www.steriking.com e-mail: steriking@wipak.com